

Cutimed® Hydro



**Hydrocolloid
dressing**

Cutimed® Hydro – focused wound care management




Cutimed® Hydro is a range of sterile, self-adhesive hydrocolloid dressings including Cutimed® Hydro L, mainly for dry to lightly exuding wounds, and Cutimed® Hydro B for light to moderately exuding wounds.

Hydrocolloids absorb exudate to form a gel, which is removed in one piece with the dressing. Cutimed® Hydro dressings promote healing by creating and controlling moist wound conditions.

Cutimed® Hydro – mode of action

The absorbent and expandable particles in the hydrocolloid layer swell as they absorb exudate and form a gel.

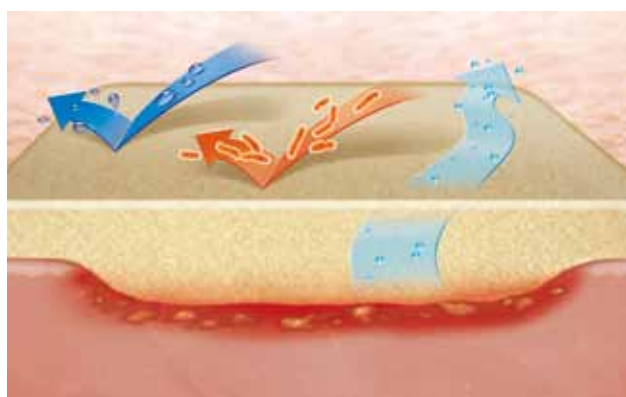

- **Facilitates moist environment**
 - promotes formation of granulation tissue
- **Forms gel in contact with exudate**
 - molds itself to structure of wound
- **Vertical absorption of exudate**
 - keeps wound edges dry and reduces risk of maceration
- **Seals the wound**
 - minimizes contamination from external fluids or bacteria
- **Product structure allows removal of the gel with the dressing**
 - minimal damage to newly forming tissue

Cutimed® Hydro – user benefits

- **Wide range of dressing shapes and sizes**
 - easy application
- **Smooth and flexible dressing**
 - conforms easily to skin and body contours
- **Intelligent dressing surface**
 - adheres securely to healthy skin, gels on wound





Cutimed® Hydro – patient comfort

- **Retains exudate in gel**
 - reduces leakage risk
 - keeps patient as active as possible
- **Waterproof dressings**
 - can be worn in the shower
- **Rapid dressing changes**
 - minimal patient discomfort
- **Free of latex and gelatin**
 - low allergy risk for patient and clinician



Cutimed® Hydro range

Hydro L is for 'Light', because the top layer of semi-permeable, low-friction polyurethane film allows close monitoring of the wound.

Hydro B is for 'Border', because the edges of the dressing are tapered for superior adhesion and minimal leakage. The top layer of the dressing is a film-coated polyurethane foam.

Indications

Cutimed® Hydro L is for dry to lightly exuding, mainly acute wounds, such as:

- Laceration
- Minor abrasions
- Minor burns (superficial and partial thickness)
- Post-op protection of closed surgical sites



Cutimed® Hydro B is for light to moderately exuding, mainly chronic wounds, such as:

- Venous ulcer
- Arterial ulcer
- Pressure ulcer
- Diabetic ulcer



NOTE: Cutimed® Hydro L can also be used as semi-occlusive secondary dressing.

Cutimed® Hydro

- Maintains a moist environment
- Stimulates granulation
- Retains exudate in gel
- Waterproof and breathable

Cutimed® Hydro L

Size (inches)	REF. No.	Dressings/Box	HCPCS
2 x 10	7263600	12	A6235
3 x 3	7263601	10	A6234
4 x 4	7263602	10	A6234
6 x 6	7263603	5	A6235

Cutimed® Hydro B

Size (inches)	REF. No.	Dressings/Box	HCPCS
3 x 3	7263500	5	A6238
4 x 4	7263501	5	A6237
6 x 6	7263502	5	A6238
Sacral 6 x 7	7263503	3	A6238

Wound Phase

-  Necrotic
-  Infected
-  Sloughy 
-  Granulation 
-  Epithelization 

Wound Depth

-  Shallow 
-  Deep
-  Shallow + deep

Exudate

-  Dry to low
-  Low to medium 
-  Medium to high



Please refer to the product label and / or package insert for full instructions on the safe use of these products.

BSN medical Inc. 5825 Carnegie Blvd. Charlotte, NC 28209-4633
Tel. (+1) 704 554 9933 Fax (+1) 704 358 4558 www.bsnmedical.com
To order toll-free: BSN medical (+1) 800 552 1157